



**SCHÖNN**  
Medizintechnik GmbH


# MEDICAL GAS ALARM MANAGEMENT SYSTEM




# MEDICAL GAS ALARM MANAGEMENT SYSTEM

## Features

- It allows comprehensive monitoring of all the medical gas alarm panels in the hospital.
- Data of the alarm panels can be displayed at a central monitoring station.
- All gas alarm panels in the hospital can work integrated with this system.
- The instantaneous pressure and alarm conditions of the gases of all gas panels in the hospital can be monitored from the central computer.
- With this system, each gas alarm level (high alarm level and low alarm level) can be changed on the gas alarm panels.
- When any alarm occurs on the gas alarm panels, the alarm can be observed both in the central panel and in this system.
- With Modbus communication system, up to 255 alarm panels can be monitored and controlled from this system.
- When any alarm occurs, the information of the panel with the fault alarm is opened as a popup.
- Alarms can be silenced remotely for each panel.
- Language select on (German-English)
- Authorized person access by password.

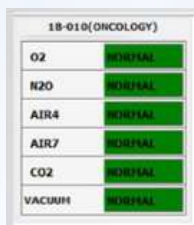

all gases will appear "normal" when there is no fault.

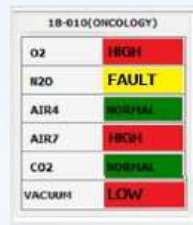

the gases will appear "high, low or fault" when there is a fault.



When any alarm occurs, the information of the panel with the fault alarm is opened as a popup.



### Mute

Alarms can be silenced remotely for each panel.

### Settings

Gas alarm level (high alarm level and low alarm level) can be changed on the gas alarm panels.

Headquarters/Factory	Warehouse
<p>📍 Andreas-Bornes Str. 46, 41179 Mönchengladbach /Germany</p> <p>☎ +49 (0) 2161 564 5075    📠 +49 (0) 2161 564 5074</p>	<p>📍 Herforder Strasse 46, 32602 Vlotho/Germany</p> <p>☎ +49 (0) 5733 963 0442    📠 +49 (0) 5733 963 0440</p>