



SCHÖNN
Medizintechnik GmbH

**PATIENT BEDHEAD
UNITS**



PATIENT BEDHEAD UNITS

PATIENT BEDHEAD UNITS



- ◆ Innovative Design
- ◆ Easy and Quick Installation
- ◆ Wide options of design, color, and dimensions
- ◆ Laminated aluminum composite panels.
- ◆ Easy to clean, antibacterial, hygienic.
- ◆ Robust support system of BHU is ideal for renovation.
- ◆ Prefabricated sockets, outlets, and assemblies with the ward's medical piping and wiring system.

The Schön bed head unit is perfectly suited for patients' rooms and other departments that require suitable lighting, electricity, telecommunications, and medical gas outlets at each bed. The unit housing is made of extruded aluminum and is powder-coated. Schön bed head units fully conform to HTM08-03, ISO EN 11197 and IEC EN 60601/1 standards, providing users with the following benefits.

DESIGN

The shape of the profile itself is a great leap forward from classic rectangular channels. The electricity and low-power outlets are also available in the most modern designs. The user or designer can choose the length of the units, select a color from the RAL chart, decide on the intensity of the lighting, and choose the number and positioning of separate outlets.

LIGHT


Shönn bedhead units can be equipped with indirect lights for room illumination, direct lighting for illuminating the bed head areas, and a night lamp. Positioning the direct lighting fixtures at a certain angle allows the light to disperse over the patient's head and provide excellent illumination for reading. Positioning of the indirect lighting efficiently disperses light over the wall above the unit. All lights are fitted with high-quality, long-life, environmentally friendly, latest-technology LED strip tubes meeting the highest standards.

SECURITY

Inside of the unit is divided into several channels. Such a layout allows a separate channel for each energy source (high voltage, low voltage, communications, medical gases), preventing any chance of interference between energy sources.

PATIENT BEDHEAD UNITS

DK 2014

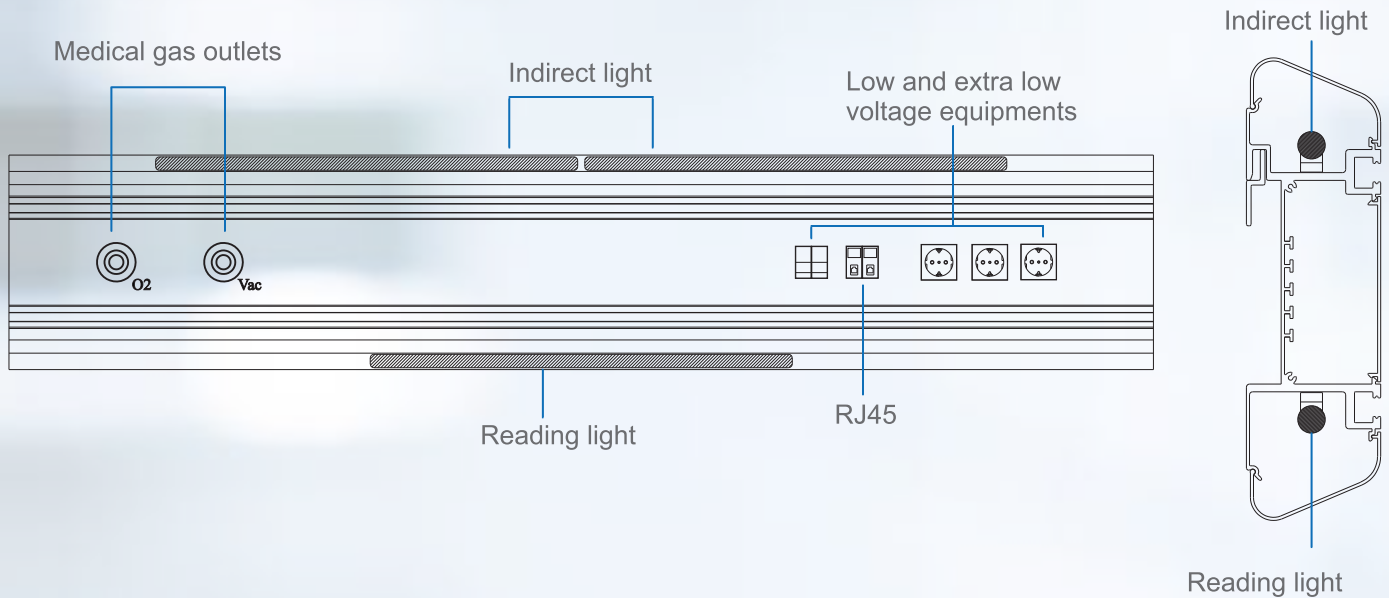



DK 2014	
Length	Product Code
100 cm	S90-750-100
120 cm	S90-750-120
150 cm	S90-750-150
160 cm	S90-750-160
180 cm	S90-750-180
200 cm	S90-750-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

SC 2015

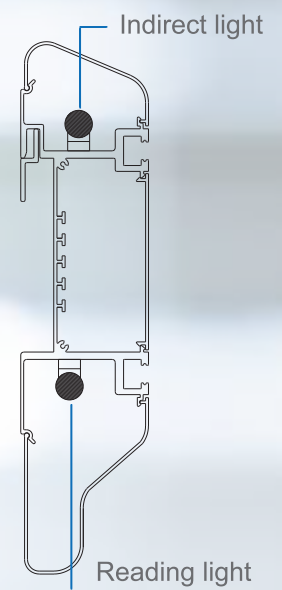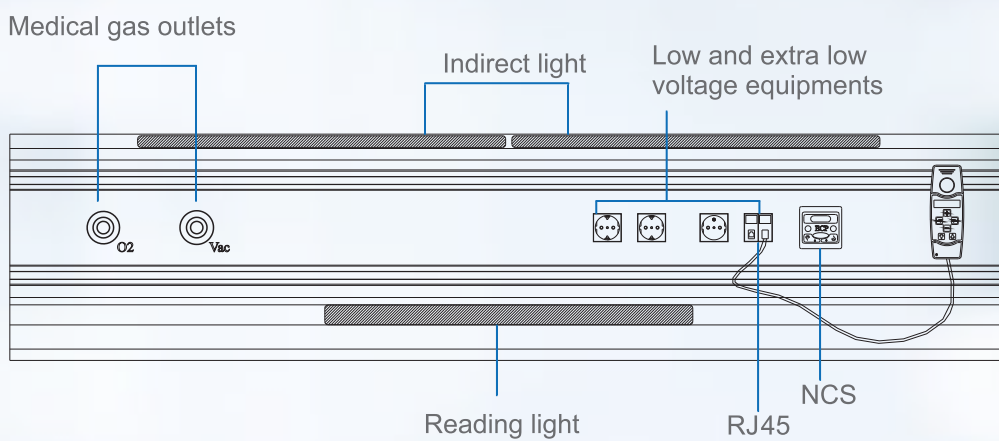

SC 2015



Length	Product Code
100 cm	S90-760-100
120 cm	S90-760-120
150 cm	S90-760-150
160 cm	S90-760-160
180 cm	S90-760-180
200 cm	S90-760-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

DK 2014-M

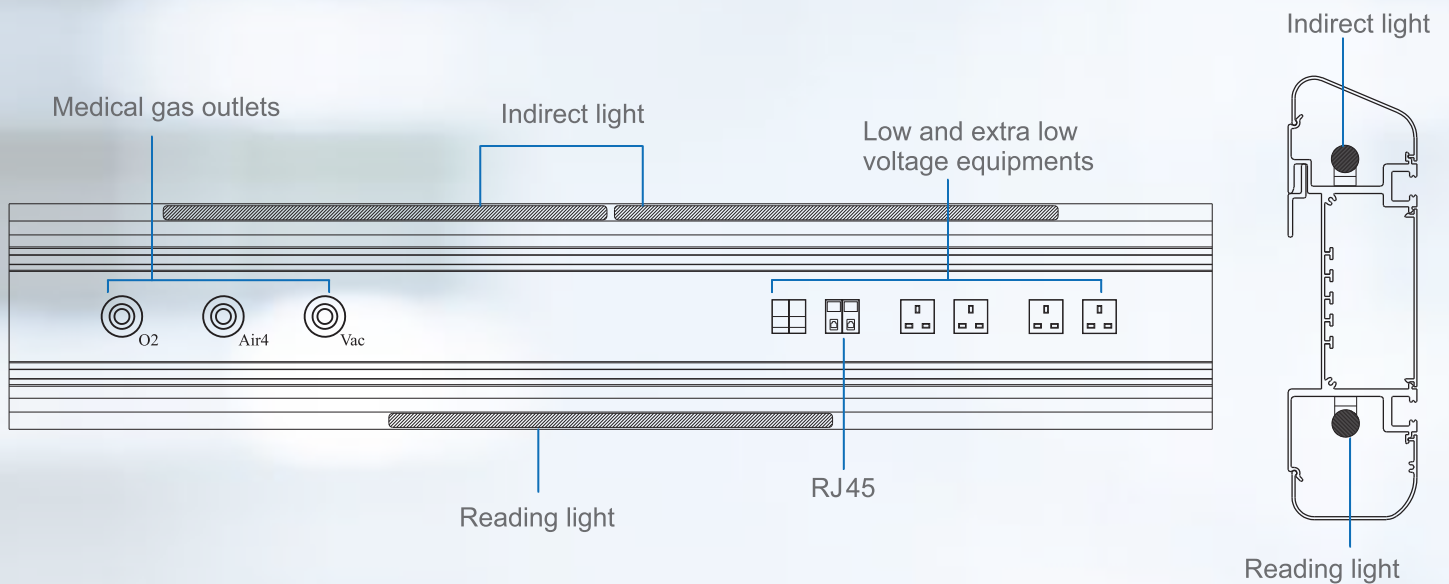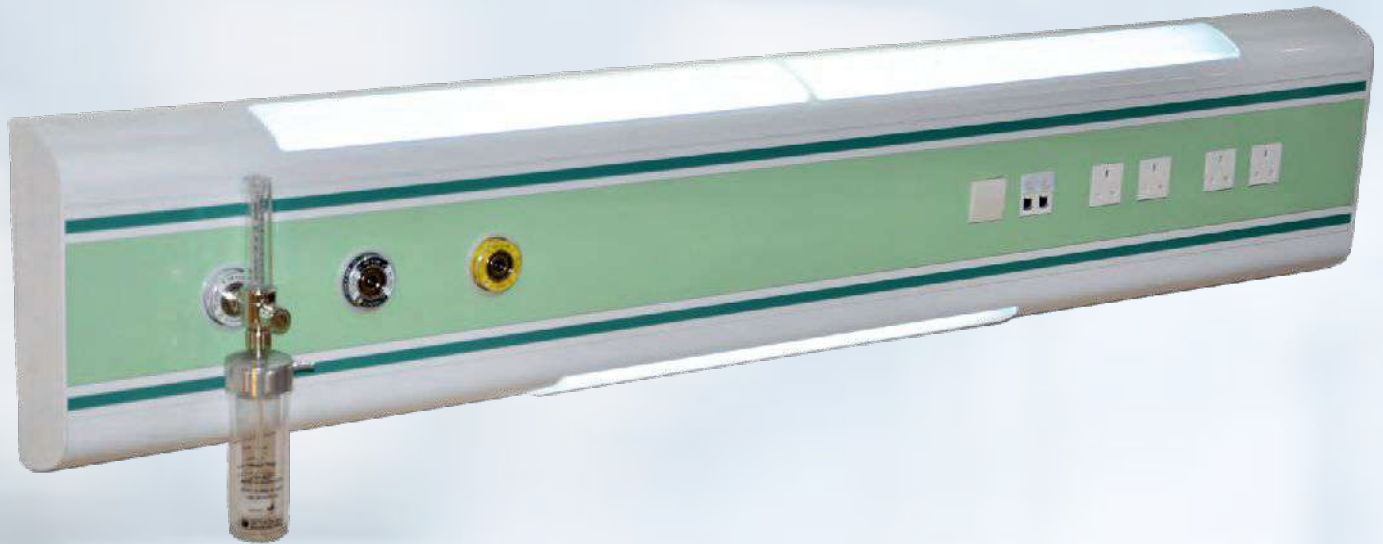

DK 2014-M



Length	Product Code
100 cm	S90-751-100
120 cm	S90-751-120
150 cm	S90-751-150
160 cm	S90-751-160
180 cm	S90-751-180
200 cm	S90-751-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

HY 2017

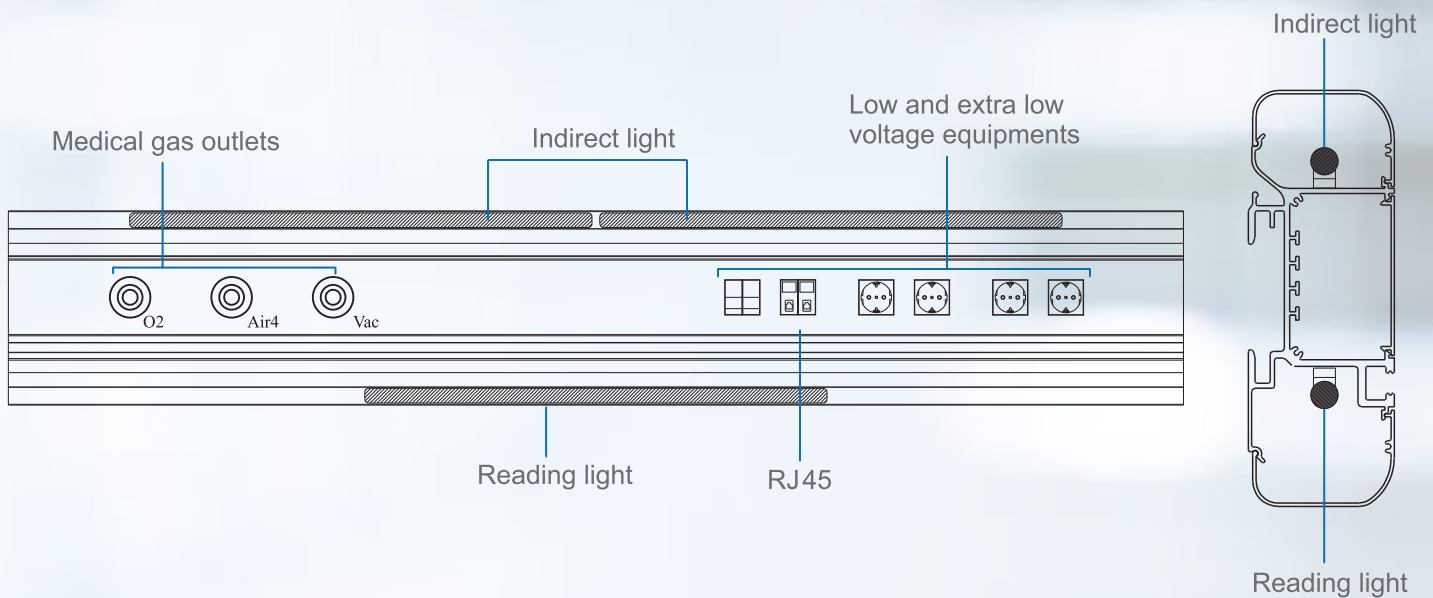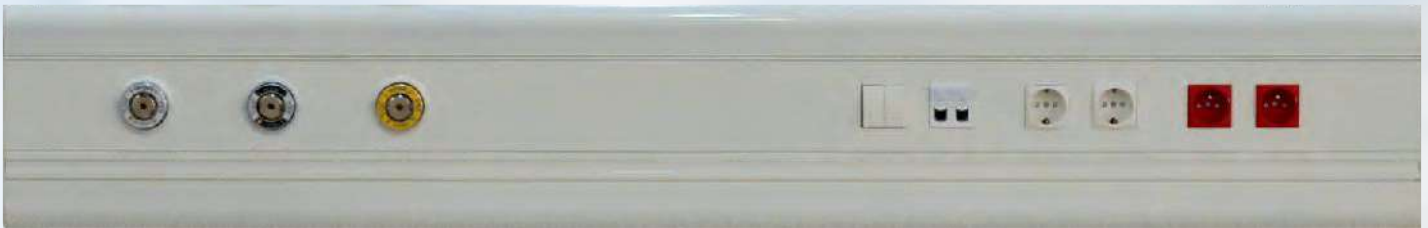

HY 2017


Length	Product Code
100 cm	S90-780-100
120 cm	S90-780-120
150 cm	S90-780-150
160 cm	S90-780-160
180 cm	S90-780-180
200 cm	S90-780-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

HY 2017-DS

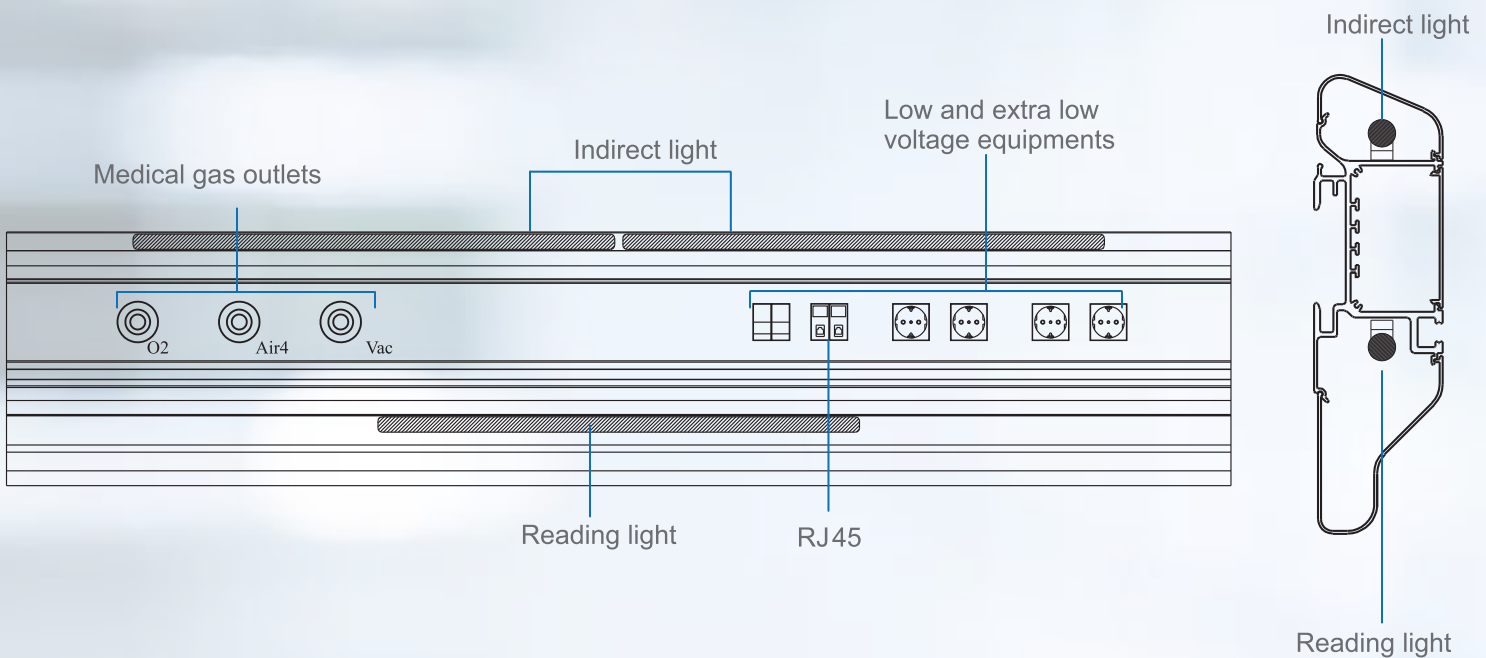

HY 2017-DS


Length	Product Code
100 cm	S90-781-100
120 cm	S90-781-120
150 cm	S90-781-150
160 cm	S90-781-160
180 cm	S90-781-180
200 cm	S90-781-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

HY 2017-D

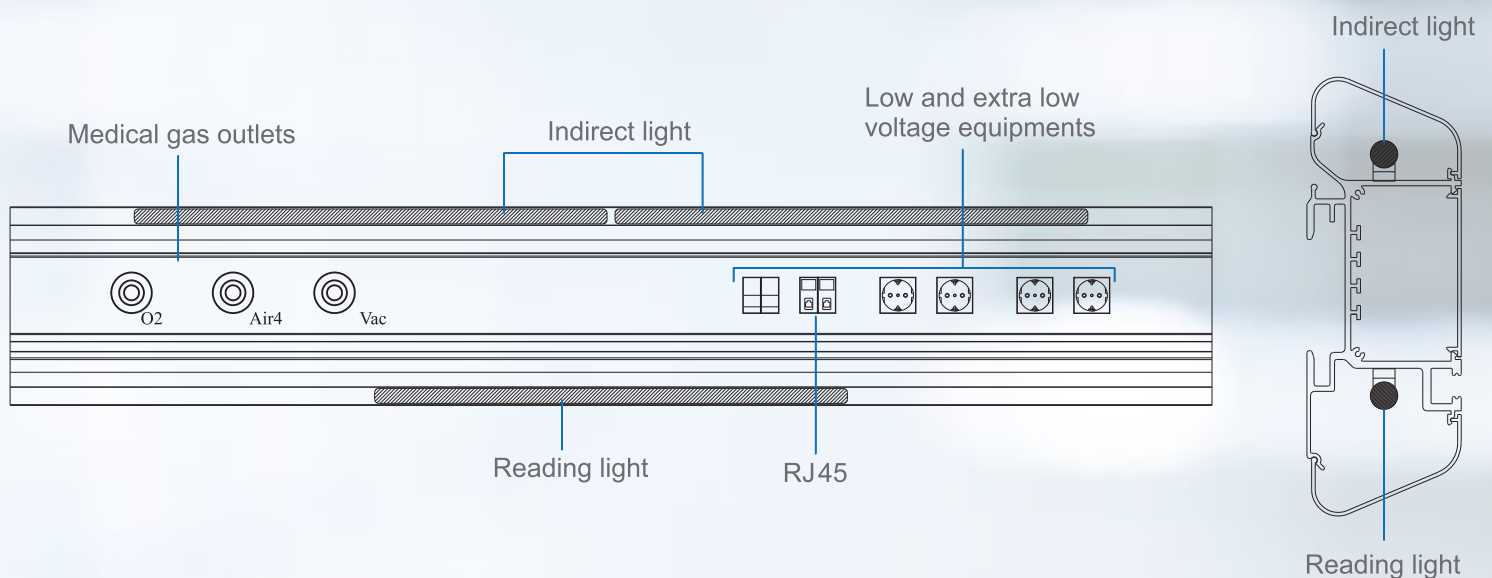

HY 2017-D



Length	Product Code
100 cm	S90-782-100
120 cm	S90-782-120
150 cm	S90-782-150
160 cm	S90-782-160
180 cm	S90-782-180
200 cm	S90-782-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.

PATIENT BEDHEAD UNITS

MD-2014

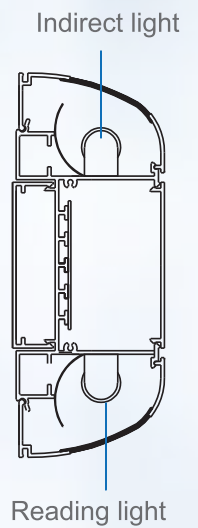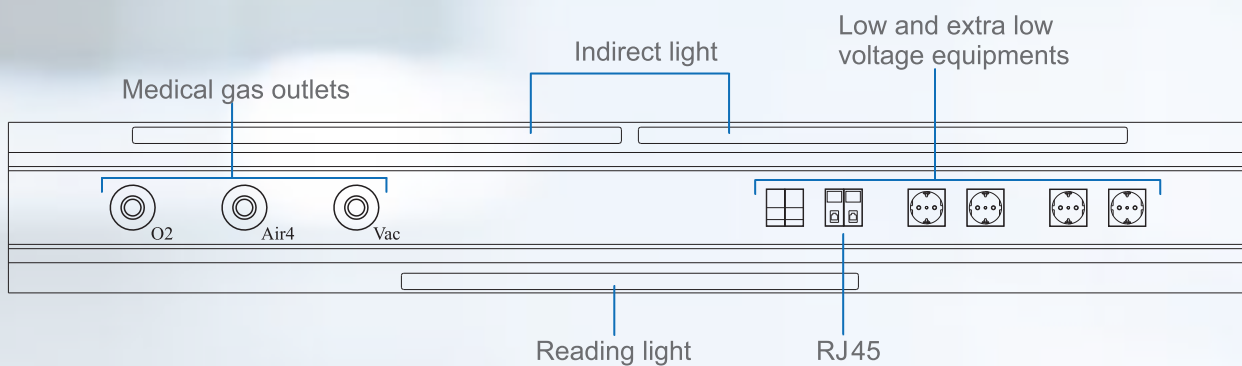

MD-2014	
Length	Product Code
100 cm	S90-320-100
120 cm	S90-320-120
150 cm	S90-320-150
160 cm	S90-320-160
180 cm	S90-320-180
200 cm	S90-320-200

	RAL 7035		RAL 7001
	RAL 3002		RAL 3014
	RAL 6019		RAL 9001
	RAL 5017		RAL 5018
	RAL 5024		Optional

Any feature of unit can be customized as per customer requests.